

# **Laserglow Product Datasheet**

# LRS-0556 DPSS Laser System

# **Laserglow Part Number: R55500XSX**

This model is listed as inactive in our product database. Stock may be limited, and availability is subject to

change without notice.

#### Similar Products: For information about the other lasers in this product family visit: http://www.laserglow.com/R55

## Ordering:

Order Online Now or Request Quote: http://www.laserglow.com/R55500XSX

## Series Specifications:

Nominal Wavelength	556 nm
Output Type	CW
Laser Source Type	DPSS



## **Overview:**

The LRS-0556 Series of Diode-Pumped Solid-State (DPSS) Lasers are ideal for applications requiring 556 nm laser light at power levels from 5 mW to over 2000 mW. This series is available in two levels of long-term output power stability and has a long operating lifetime at an aggressively competitive cost.

These lasers are commonly used for Raman Spectroscopy and other applications requiring a narrow linewidth and good long-term output power stability. The driver is available as a complete FDA-compliant system or as an O.E.M. component with significantly reduced dimensions.

Laserglow products are currently being used by some of the World's top universities and other prominent research facilities.

#### **Key Features:**

- · Air cooled no need for water cooling or external chiller
- Lightweight, compact design
- Efficient DPSS technology runs on standard AC power (85 264 V, 47 63 Hz)
- >10,000 hours continuous maintenance-free operating life
- TTL and Analog modulation (input via BNC connector) lab-spec models only
- Adjustable output power (via lockable dial) lab-spec models only
- LED display showing LD current, laser cavity temperature lab-spec models only
- FDA CDRH Compliant Class IIIb / Class IV enclosure
- · 48-hour replacement coverage available for an additional fee on specific models

- Laser Head
  Driver/Power Supply
  Power Cable
  BNC Connector (LabSpec models only)
  Keys, Safety Interlock
  Hard-shell Carrying Case

## **Specifications:**

This spec sheet has been generated specifically for part number R55500XSX, per your request, and data for the entire series is also displayed for your reference. The specs which are specific to R55500XSX have been highlighted below in red + bold.

Output Power (mW)	>30, >50, >80, >100, >200	>300, <b>&gt;500</b> , >800, >1000	>1500, >2000, >2500
Output Power Stability (%RMS/4h)	<1, <3, <5	<1, <3, <5, <mark>&lt;10</mark>	<3, <5
Central Wavelength (nm)	556.3	556.3	556.3
Wavelength Tolerance (+/- nm)	1	1	1
Divergence (mrad, full angle)	<1.5	<1.5	<2
Beam Dimensions (mm, 1/e²)	2	4	4
Warm-up Time (minutes)	10	10	10
Spectral Linewidth (nm)	<0.15	<0.15	<0.15
M <sup>2</sup>	<1.2	<3	<6
Polarization Ratio	>100	>100	
Beam Pointing Stability (mrad)	<0.05		
Operating Temperature Range (°C)	10 to 35	10 to 35	10 to 35
Max. Analog Modulation Freq. (Hz)	500	500	500
Max. TTL Modulation Freq. (Hz)	500	500	500
Modulation Input Signal	0-5 VDC	0-5 VDC	0-5 VDC
Total Power Consumption (W)	30, 70		
Max. Power Input Duty Cycle	100%	100%	100%
Standard Warranty (months)	12	12	12
MTTF (operational hours)	10000	10000	10000
Weight of Product or Laser Head (kg)	2	2.6	6.1
Beam Height from Base Plate (mm)	27.4	68.2	93.5
Dimensions of Product or Laser Head (mm)	197 (l) x 70 (w) x 50 (h)	240 (l) x 99 (w) x 94 (h)	346 (I) x 140 (w) x 125 (h)

CW: All specifications are based on performance at full output power and after the specified warmup period. Output characteristics may change if the laser is run at a different power level.

Q-Switched: Specifications are based on the laser pulsing at the specified design frequency. Output characteristics may change if the laser is run at a

different frequency.

## **Power Supply Options:**

These lasers are available with several different power supply options. The model that you have selected is highlighted below, and any other options are shown for easy reference.

	Power Supply Type:	FO	FN	FW
FDA-Compliant LabSpec	Input Power	85v to 264v	85v to 264v	85v to 264v
3 mg 2-	Power Supply Weight (kg)	2.6	2.6	5.2
and the second sec	Dimensions (mm)	268 (l) x 145 (w) x 106 (h)	268 (l) x 145 (w) x 106 (h)	307 (l) x 168 (w) x 123 (h)


	Power Supply Type:	SN
FDA-Compliant Standard	Input Power	85v to 264v
	Power Supply Weight (kg)	2.3
	Dimensions (mm)	238 (l) x 146 (w) x 102 (h)

\*Power supply may not be exactly as shown, see dimensional drawings on next 2 pages. \*Dimensions for fiber-integrated (I\_) include laser head packaged inside.


## **Regulatory Classification:**

The model you have selected (R55500XSX) requires the following safety label(s):



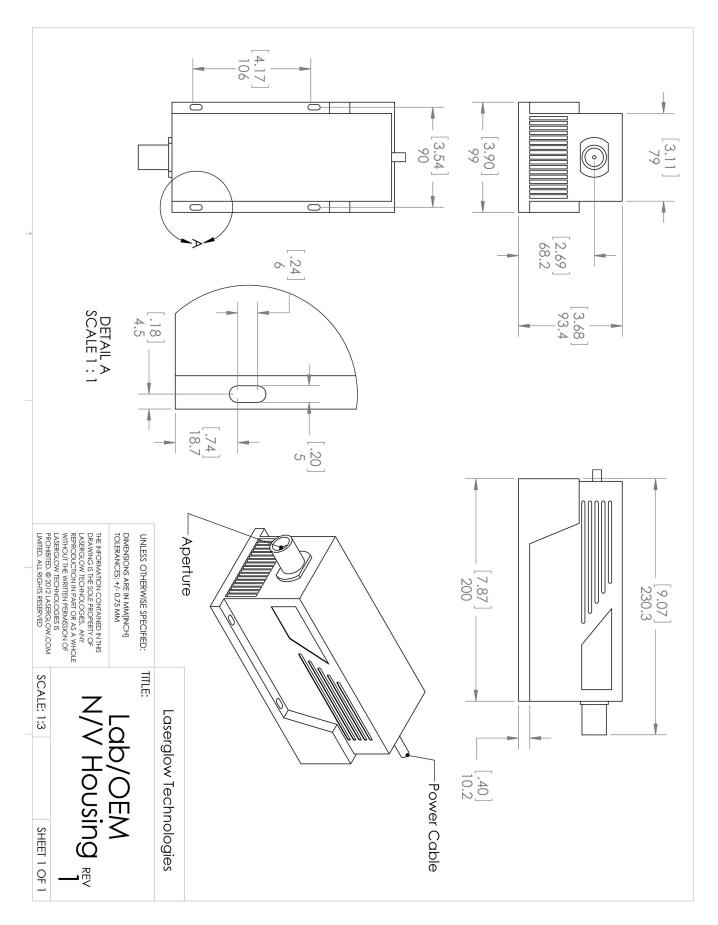



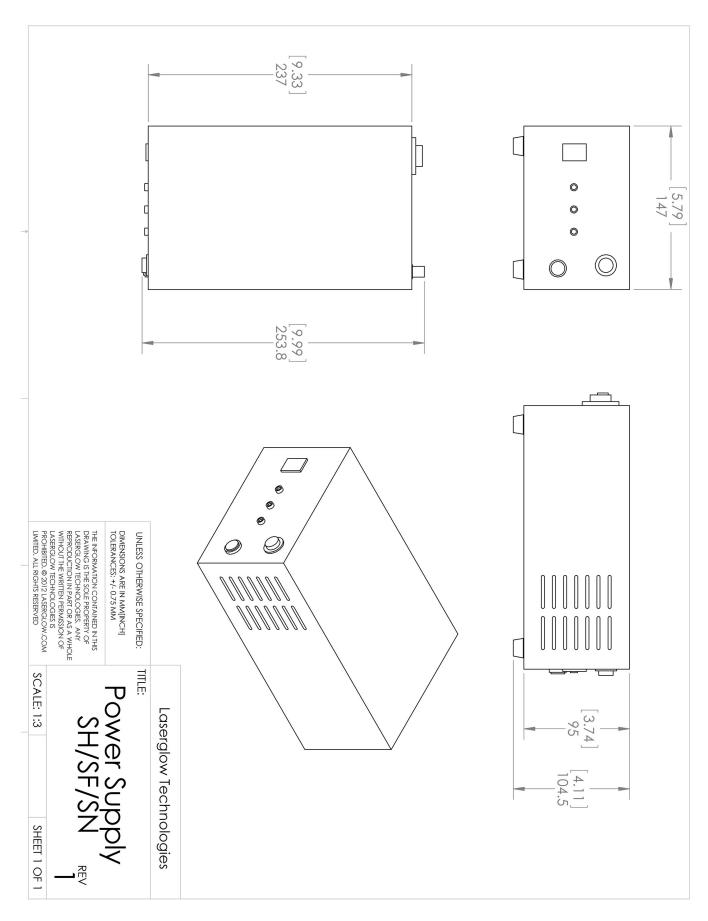

## Accessories:

The most popular accessories for model R55500XSX are shown below. For additional details regarding these or other accessories please see our website or contact us directly.

Part Number	Description	
AGF5565XX	LSG-556-NF-5 Fit-Over Safety Goggles 556nm Output: OD 5+ at 190-579 nm CE Certified Full Details: www.laserglow.com/AGF	
ACFVISHXA	FC/PC Fiber Coupler/Collimator for visible spectrum wavelengths (400 to 700 nm) (installed and aligned) 11mm diameter input lens Full Details: www.laserglow.com/ACF	
ACSVISHXA	SMA-905 Fiber Coupler/Collimator for visible spectrum wavelengths (400 to 700 nm) (installed and aligned) 11mm diameter input lens Full Details: www.laserglow.com/ACS	
AFS2002XX	Armored Fiber With SMA 905 Connectors 200um Core Multimode 2 m length Full Details: <u>www.laserglow.com/AFS</u>	
AFF2002XX	Armored Fiber With FC/PC Connectors 200um Core Multimode 2m length Full Details: <u>www.laserglow.com/AFF</u>	
ТВК	Complete optics kits with breadboard mounting hardware. External modulators, variable attenuators, free-space fiber launch systems Full Details: <u>www.laserglow.com/TBK</u>	
ACALBNWXX	Carrying Case-104 Holds Lab/OEM N or W size Standard/Labspec laser Full Details: <u>www.laserglow.com/ACA</u>	Included With Laser

## FOR MORE INFORMATION PLEASE CONTACT:

LASERGLOW TECHNOLOGIES 99 Ingram Dr. Unit B, North York, ON, Canada M6M2L7 Tel. (416) 729-7976 Fax (716) 322-3510 sales@laserglow.com www.laserglow.com

E&OE: Data included in this sheet may be subject to change without notice. Please confirm critical specifications with our staff prior to ordering.